




MEDEQUIP
HOME CARE SERVICES

UNIVERSITY OF MICHIGAN
HEALTH SYSTEM


Oxygen Concentrator Tubing Assembly, with Humidification

How do I connect my nasal cannula, tubing and humidifier to my concentrator?

Step 1: Connect your nasal cannula to the tubing swivel connector

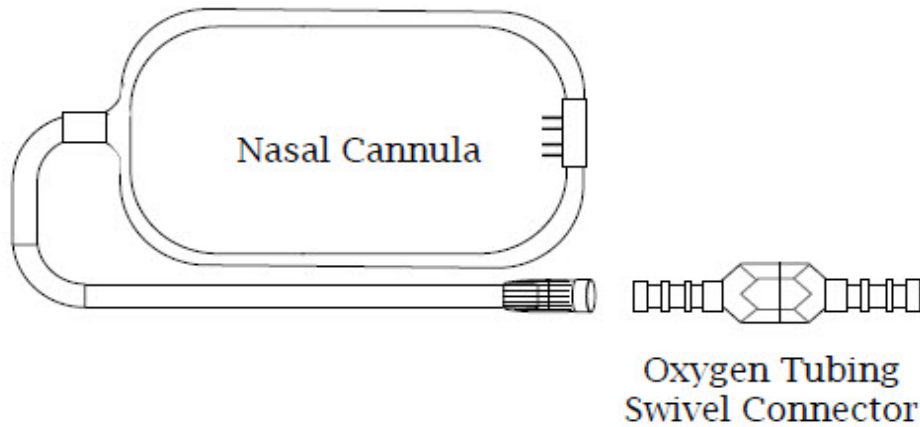


Step 2: Connect your swivel connector to your tubing

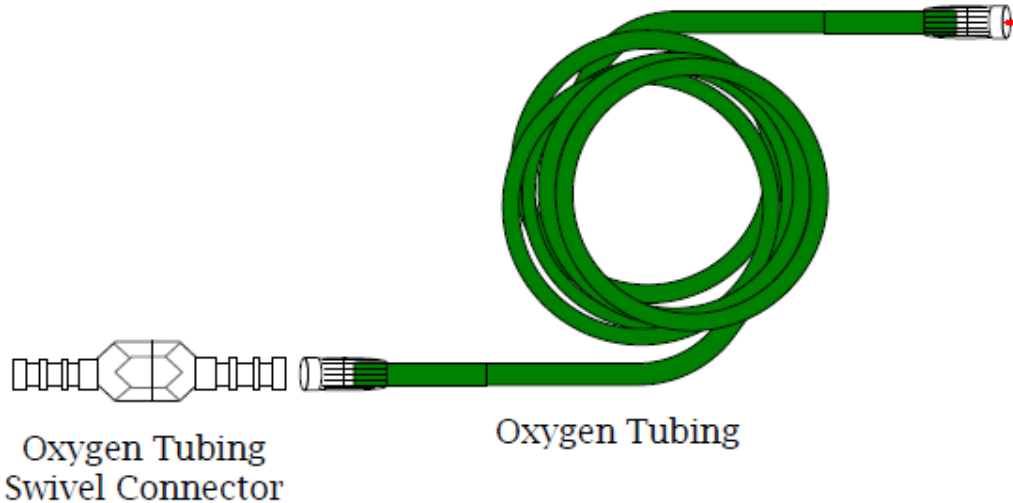


Home Care Services, MedEQUIP
(800) 530-0714

Step 3: Connect your oxygen tubing to the humidifier

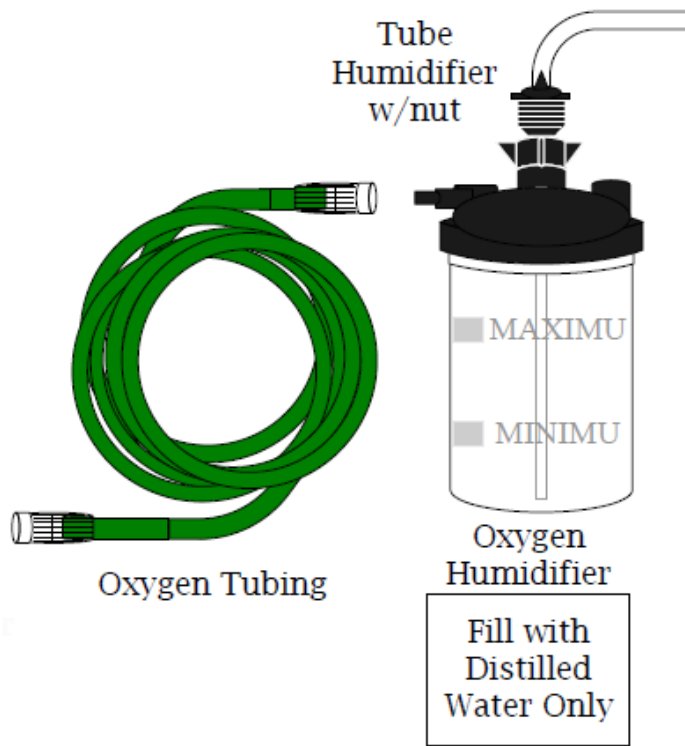


Step 4: Connect the humidifier with the nut to the concentrator unit

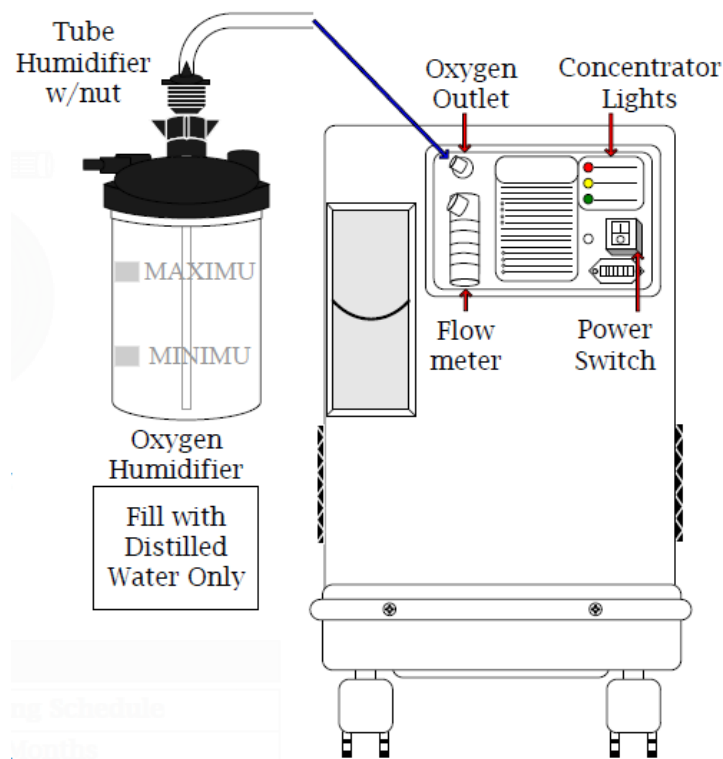

Home Care Services, MedEQUIP
(800) 530-0714

Oxygen Concentrator Assembly with Humidification - Adult

Step 5: Total unit connection

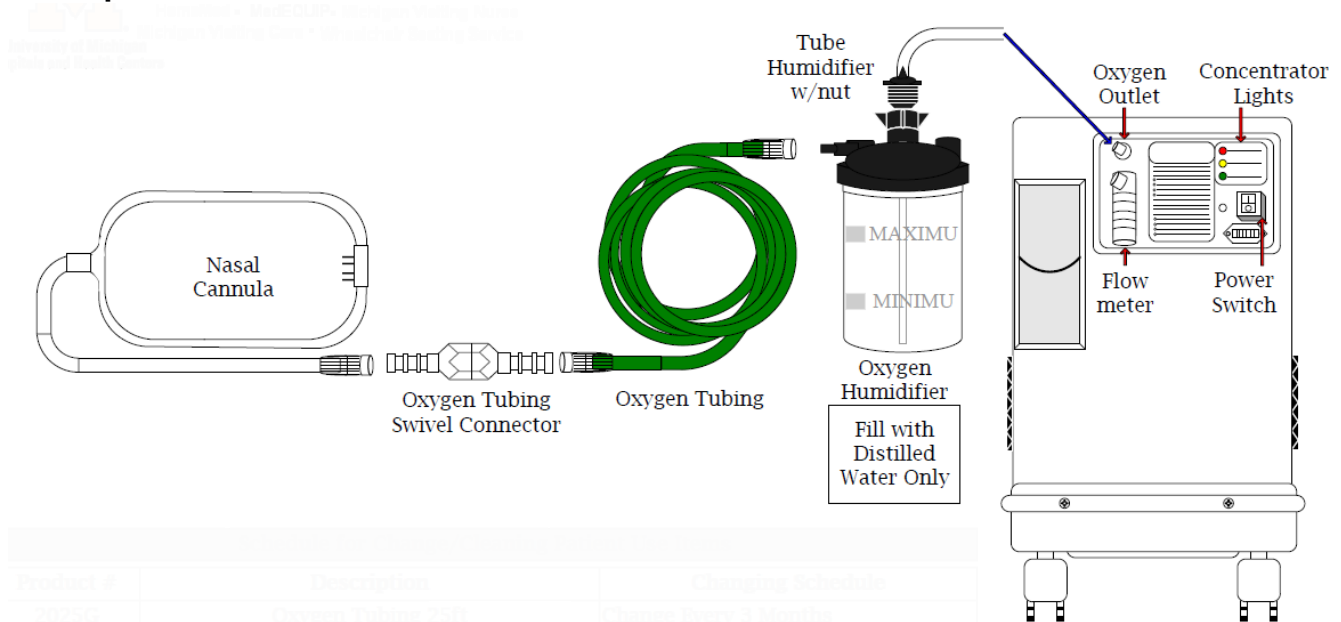


How often should I change my supplies?

Schedule for Change/Cleaning Patient Use Items

Product #	Description	Changing Schedule
2025G	Oxygen Tubing 25ft	Change Every 3 Months
2050G	Oxygen Tubing 50ft	Change Every 3 Months
1600/1601	Nasal Cannula	Change Every 2 Weeks
1220	Oxygen Tubing Swivel Connector	Change Every 3 Months
Wipe and clean nasal cannula with damp cloth twice daily		
7600	Oxygen Humidifier	Change Every 2 Weeks
	Tube Humidifier w/Nut	Only Replace if Damaged
Daily Wash with Warm Water using Soap. Weekly, Soak in 1 part White Vinegar to 3 Parts Water for 30 minutes Rinse thoroughly and allow Parts to Air Dry.		

Copyright © 2015 The Regents of the University of Michigan, Illustrated By Rachel Stout

What indicator lights do I need to know on my concentrator?

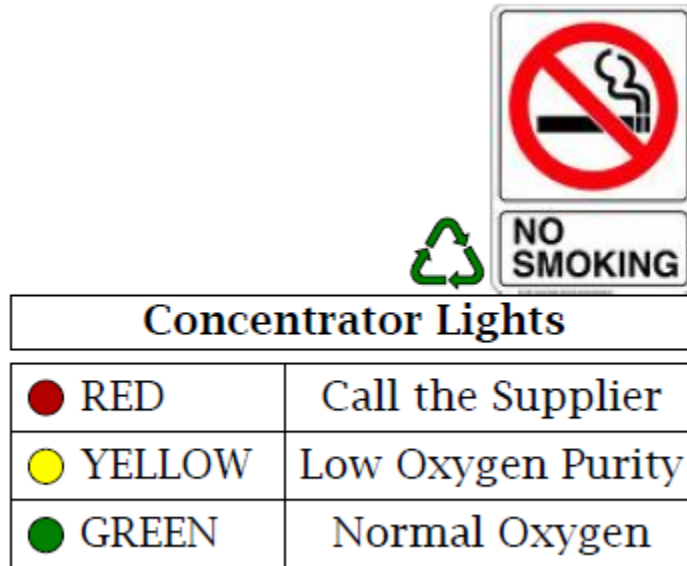

Disclaimer: This document contains information and/or instructional materials developed by the University of Michigan Health System (UMHS) for the typical patient with your condition.

It may include links to online content that was not created by UMHS and for which UMHS does not assume responsibility. It does not replace medical advice from your health care provider because your experience may differ from that of the typical patient. Talk to your health care provider if you have any questions about this document, your condition or your treatment plan.

Author: Rachel Stout, Garret Fugitt, RRT, LRT, Brad Hayden

Reviewers: Sarah Middlemas, MPH

Patient Education by [University of Michigan Health System](#) is licensed under a [Creative Commons Attribution-NonCommercial-ShareAlike 3.0 Unported License](#). Last Revised 02/2016