


Fiasp® FlexTouch® (insulin aspart injection) 100 units/mL


Fiasp FlexTouch Insulin Pen parts:

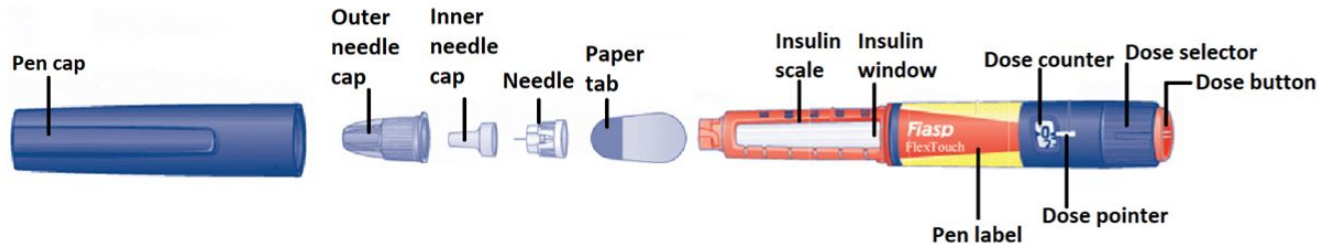


© 2020 Novo Nordisk Canada Inc. All rights reserved

Where to inject:

- Rotate injections between injection spots #1-4 in the diagram to the right:
 1. **Stomach**, except for the 2-inches around the belly button
 2. Top and outside of **thighs**, except in the bony area close to the knees
 3. Back of the **upper arms** (might work best if someone else is injecting)
 4. Top and outside of the **buttocks**
- Insulin is injected into the layer of fatty tissue under the skin, not into the muscle
- Do **not** inject into the same spot every time
- Do **not** inject into moles or scars

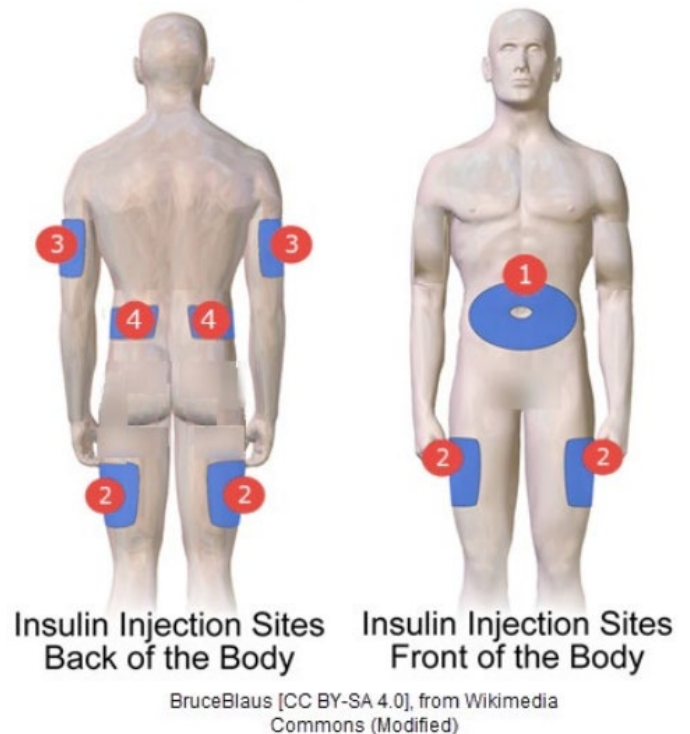

Step 1 Check The insulin

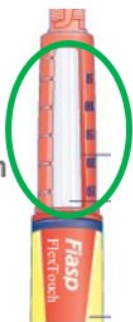

Prepare to inject:

1. If your Fiasp FlexTouch pen is in the refrigerator, take it out at least 1-2 hours before using
2. Check the label to make sure you are using the right insulin pen
3. Check the expiration date
4. Remove the pen cap to check the insulin
 - Do **not** use insulin that is cloudy or not clear, or if it has any specks



Step 2
Attach
The needle

5. Wipe the pen tip with an alcohol pad
6. Remove the seal on the needle cap
 - Do **not** share needles or insulin pens with anyone else
 - Do **not** re-use needles
7. Twist or push (based on the needle type) the needle straight onto the pen tip

Step 3
Perform
An air shot




For each injection:

1. Select a dose of 2 units
2. Take off the outer needle cap (save it) and inner needle cap (throw it away)
3. With the pen pointing up, tap the insulin to move the air bubbles to the top



4. Press the button all the way in and make sure insulin comes out of the needle (this is called an air shot)
 - Repeat up to two more times with the same needle if needed
 - If insulin does not come out after three times, change the needle and try again
 - If insulin still does not come out after changing the needle, the pen may be broken

Step 4
Select
The dose

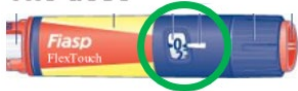

5. Check that the dose counter shows “0” after the safety test
6. Turn the dose counter to the number of Fiasp FlexTouch units that equals your dose
 - Do **not** push the button while selecting your dose

Step 5
Inject
The dose



7. Clean the injection site with an alcohol swab and wait for it to dry
8. Put the needle into the skin all the way
9. Press and hold the button to give the dose
10. Keep the button pressed and slowly count to 10 before taking the needle out of the skin
 - Do **not** rub the spot where you injected

Step 6
Remove
And discard
needle



11. Put the outer needle cap back on the needle and use it to remove the needle from the pen
 - Remove the needle from the pen after each dose
12. Place used needles in red sharps container or a hard plastic container (ex. an empty laundry detergent bottle)
 - Store the Fiasp FlexTouch pen without a needle attached and with the pen cap on

Where to store Fiasp FlexTouch:

- Keep unused pens in the refrigerator
- Keep open or in-use pens at room temperature for a maximum of 28 days
- After 28 days, Fiasp FlexTouch pens should be thrown away in the trash - even if there is leftover insulin
- Do **not** keep Fiasp FlexTouch where there is direct heat or direct light
- If Fiasp FlexTouch is frozen or overheated, throw it away

Top five tips for using your insulin:

1. Use Fiasp FlexTouch up to 2 minutes before, or 20 minutes after starting a meal or snack, as directed
2. Rotate your injections between the four injection spots listed above
3. Always use a new needle for each injection and clean the injection site before each dose
4. Slowly count to 10 after pressing the dose button
5. If insulin leaks after the injection, do **not** re-dose; be sure to keep the needle in the skin for at least 10 seconds next time

If you would like to watch a video on how to inject insulin using your Fiasp FlexTouch, you can visit:

<https://www.rapidactinginsulin.com/fiasp/using-fiasp/videos-and-downloads.html>



Disclaimer: This document contains information and/or instructional materials developed by Michigan Medicine for the typical patient with your condition. It may include links to online content that was not created by Michigan Medicine and for which Michigan Medicine does not assume responsibility. It does not replace medical advice from your health care provider because your experience may differ from that of the typical patient. Talk to your health care provider if you have any questions about this document, your condition or your treatment plan.

Edited by: Karelyn Munro, BA

Patient Education by [Michigan Medicine](#) is licensed under a [Creative Commons Attribution-NonCommercial-ShareAlike 4.0 International Public License](#). Last Revised 3/2020